



## Hilsonic Ultrasonic Cleaning Equipment

Units 1-4, Cashel Road, Wirral,  
Merseyside CH41 1DY

Tel: +44 (0)151 639 6020  
Fax: +44 (0) 151 334 7407

[www.hilsonic.co.uk](http://www.hilsonic.co.uk)



### EXAMPLE APPLICATIONS

Dental Clinics  
Medical Rooms  
Optical Laboratories  
Chemical Laboratories  
Electrical Component Makers  
Watch / Metal Part Makers  
Shooting Clubs  
Special Education Institutes  
Golf Clubs  
Jewellery Makers  
Mobile Phone Service Shops  
Printing Shops  
Automotive Repair Shops

### SPECIFICATION

**Model Name:**  
8061 TOUCH

**Timer Settings:**  
5 to 30 minute full range timer

**Tank Size:**  
3.0L (MAX. fluid level 2.5L)

**Ultrasonic Frequency:**  
35,000Hz

**Power Consumption:**  
170W

**Power Supply:**  
(AC 220~240V 50/60Hz)

**Mains Fuse:**  
3A

**Unit Dimensions (LxWxH):**  
455x275x210mm

**Weight (Net):**  
5.1kg (11.3lbs)

**Weight (Gross):**  
6.1kg (13.4lbs)

# 8061 TOUCH

## PROFESSIONAL ULTRASONIC CLEANER

### CLEANING MADE SIMPLE for professional use

- Large 3L stainless steel tank
- Super-sized industrial grade ultrasonic transducer (80W)
- Ceramic heater (65°C)
- Touch screen control panel
- Full CE and RoHS approval
- Two cooling fans
- Drainage valve
- Quiet operation
- Engineering grade plastic

Ultrasonic cleaners function by producing sound waves that are transmitted into the tank and cleaning solution. These waves create millions of microscopic bubbles, which collapse or 'implode', releasing large amounts of energy, which scrub the surface clean. This process is called 'Cavitation'. These tiny bubbles create incredible cleaning results not only to the surface, but by penetrating the finer detail and recesses in surgical instruments and small items to ensure a professional clean prior to disinfection.

The **8061 TOUCH** professional cleaner has a super-sized ultrasonic transducer that is 30% more powerful than our standard transducers, to enable uniform distribution of ultrasonic waves. The touch screen control panel of the **8061 TOUCH** has a simple to use 5 to 30 minute full range digital timer. Set your required cleaning time and temperature and let the cleaner do the hard work for you. The **8061 TOUCH** has an easy to use drainage valve, making it simple to empty the dirty water from the tank.

The **8061 TOUCH** is ideal for the medical world due to new capacitive sensing touch screen technology which ensures easy usage and maintains a clean hygienic control panel surface in a clinical environment.

It is also possible to suspend a glass beaker within the tank providing in-direct cleaning for items that require isolated cleaning.




#### Included Accessories

Draining Hose




Plastic Basket

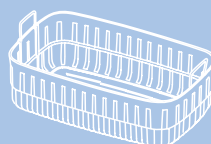

#### Optional Accessories

Stainless steel rack and tray

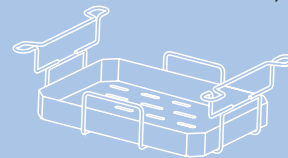

Beaker and beaker holder

